


Ordering Information
Web www.encorbio.com
Email admin@encorbio.com
Phone 352-372-7022
Fax 352-372-7066

HGNC Name: GAPDH
UniProt: P04406
RRID: AB_2572289
Immunogen: Full length recombinant human GAPDH expressed in and purified from *E. coli*
Format: Antibody is supplied as an aliquot of serum plus 5mM NaCl
Storage: Store at 4°C for short term, for longer term storage store at -20°C. Avoid freeze / thaw cycles.
Recommended dilutions:
WB: 1:30,000 IF/ICC and IHC: 1:2,000


References:

1. Morgenegg G, et al. Glyceraldehyde-3-phosphate dehydrogenase is a nonhistone protein and a possible activator of transcription in neurons. *J. Neurochem.* 47:54-62 (1986).
2. Schulze H, et al. Rat brain glyceraldehyde-3-phosphate dehydrogenase interacts with the recombinant cytoplasmic domain of Alzheimer's beta-amyloid precursor protein. *J Neurochem.* 60:1915-22 (1993).
3. Burke JR, et al. Huntingtin and DRPLA proteins selectively interact with the enzyme GAPDH. *Nature Med.* 2:347-50 (1996).
4. Dastoor Z, Dreyer J-L. Potential role of nuclear translocation of glyceraldehyde-3-phosphate dehydrogenase in apoptosis and oxidative stress. *J. Cell Sci.* 114:1643-53 (2001).
5. Hara MR, et al. S-nitrosylated GAPDH initiates apoptotic cell death by nuclear translocation following Siah1 binding. *Nat. Cell Biol.* 7:665-74 (2005).
6. Yang J-S, et al. GAPDH inhibits intracellular pathways during starvation for cellular energy homeostasis. *Nature* 561:263-67 (2018).

Applications	Host	Isotype	Molecular Wt.	Species Cross-Reactivity
WB, IF/ICC, IHC	Rabbit		36kDa	Hu, Rt, Ms, Co, Pi, Ho, Ck

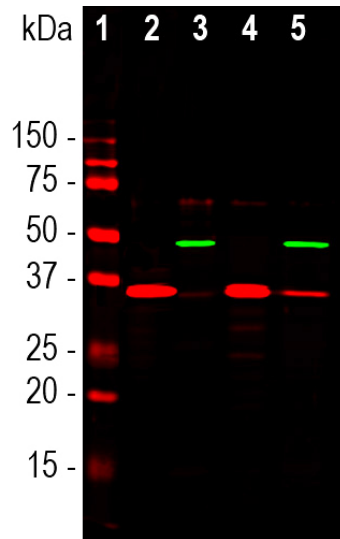

Western blot analysis of different cell cytosolic or nuclear enriched fractions, using rabbit pAb to GAPDH, RPCA-GAPDH, dilution 1:20,000, in red: [1] protein standard, [2] NIH-3T3 cytosolic, [3] NIH-3T3 nuclear, [4] HeLa cytosolic, and [5] HeLa nuclear fractions. Strong band at 37kDa corresponds to GAPDH protein, mainly detected in the cytosolic fractions. The same blot was simultaneously probed with mouse mAb to the nuclear RNA binding protein SF3B4, [MCA-3A1](#), dilution 1:1,000, in green. In contrast to GAPDH, SF3B4 is exclusively expressed in the nuclear fraction.

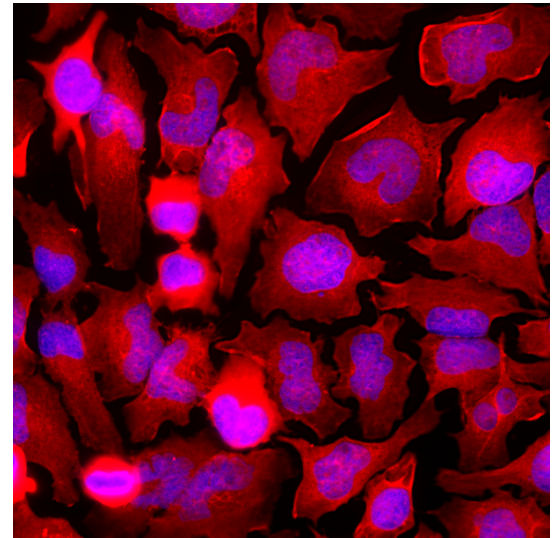

Immunofluorescent analysis of HeLa cells stained with rabbit pAb to GAPDH, RPCA-GAPDH, dilution 1:2,000 in red. Blue is Hoechst staining of nuclear DNA. The RPCA-GAPDH antibody produces diffuse cytoplasmic staining of cells.

Background:

Glyceraldehyde 3-phosphate dehydrogenase (GAPDH) is a metabolic enzyme responsible for catalyzing one step in the glycolytic pathway, the reversible oxidative phosphorylation of glyceraldehyde 3-phosphate to produce glyceraldehyde 1-3-bisphosphate. Because GAPDH protein is expressed in large amounts and is required at all times for an important "house keeping" function, levels of GAPDH mRNA are often measured and used as standards in studies of the expression of other mRNAs. Similarly specific antibodies to GAPDH are used to measure GAPDH expression as a protein standard in western blotting experiments. This antibody can be used as a loading control for such experiments, allowing comparison between the expression level of this protein and others in a cell or tissue. Apart from a role in glycolysis, GAPDH may have other roles, such as in the activation of transcription (1). GAPDH is also reported to bind to a variety of other proteins, including the amyloid precursor protein, mutations in which cause some forms of Alzheimer's disease, and the polyglutamine tracts of Huntingtin, the protein product aberrant forms of which are causative of Huntington's disease (2,3). The protein may also have a role in the regulation of apoptosis, and interestingly migrates from the cytoplasm into the nucleus when cells become apoptotic (4). The control of this is mediated by NO mediated S-nitrosylation of GAPDH in the cytoplasm which then binds to and stabilizes the Siah1 E3 ubiquitin ligase which enters the nucleus along with nitrosylated GAPDH. Nuclear Siah1 then initiates specific protein degradation events leading to the activation of a nuclear cell death pathway (5). Recently GAPDH was shown to be a substrate for the 5'-AMP dependent protein kinase (AMPK) which results in the redistribution of cellular GAPDH to cytosolic membranes and the inhibition of intracellular transport (6). AMPK is activated as AMP levels increase in response to starvation, and the pathway including GAPDH and transport inhibition apparently conserves cellular energy.

The immunogen used to raise RPCA-GAPDH was extensively purified pig GAPDH. The antibody has been widely used as a western blotting standard by many labs for many years and is known to detect GAPDH in a variety of mammalian species including human, rat and mouse. It produces a single clean band on western blots of typical cell and tissue extracts We also market a widely used mouse monoclonal antibody to GAPDH, [MCA-1D4](#) with similar properties to RPCA-GAPDH.

FOR RESEARCH USE ONLY. NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE.

Abbreviation Key:

mAb—Monoclonal Antibody **pAb**—Polyclonal Antibody **WB**—Western Blot **IF**—Immunofluorescence **ICC**—Immunocytochemistry
IHC—Immunohistochemistry **E**—ELISA **Hu**—Human **Mo**—Monkey **Do**—Dog **Rt**—Rat **Ms**—Mouse **Co**—Cow **Pi**—Pig **Ho**—Horse **Ch**—Chicken
Dr—*D. rerio* **Dm**—*D. melanogaster* **Sm**—*S. mutans* **Ce**—*C. elegans* **Sc**—*S. cerevisiae* **Sa**—*S. aureus* **Ec**—*E. coli*.

FOR RESEARCH USE ONLY. NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE.

Abbreviation Key:

mAb—Monoclonal Antibody pAb—Polyclonal Antibody WB—Western Blot IF—Immunofluorescence ICC—Immunocytochemistry IHC—Immunohistochemistry E—ELISA Hu—Human Mo—Monkey Do—Dog Rt—Rat Ms—Mouse Co—Cow Pi—Pig Ho—Horse Ch—Chicken Dr—D. rerio Dm—D. melanogaster Sm—S. mutans Ce—C. elegans Sc—S. cerevisiae Sa—S. aureus Ec—E. coli.